



**AESKULISA<sup>®</sup>**  
THE DIAGNOSTIC TOOL THAT WORKS

**INSTRUCTION  
MANUAL**

**AESKULISA<sup>®</sup> Histone-H2B**

*Ref 3153*







Product Ref.	3153
Product Desc.	Histone-H2B
Manual Rev. No.	004 : 2024-12-13

# Instruction Manual

## Table of Contents

---

1	Intended Use .....	1
2	Clinical Application and Principle of the Assay .....	1
3	Kit Contents .....	2
4	Storage and Shelf Life .....	2
5	Precautions of Use.....	3
6	Sample Collection, Handling and Storage.....	4
7	Assay Procedure .....	4
8	Quantitative and Qualitative Interpretation .....	7
9	Technical Data.....	8
10	Performance Data.....	8
11	Literature .....	9
12	Regulatory Symbols.....	10



Aesku.Diagnostics GmbH & Co. KG  
Mikroforum Ring 2  
55234 Wendelsheim  
Germany  
Phone: +49 6734 9622-0  
Fax: +49 6734 9622-2222  
Website: [www.aesku.com](http://www.aesku.com)  
Mail: [info@aesku.com](mailto:info@aesku.com)

## 1 Intended Use

**AESKULISA® Histone-H2B** is a solid phase enzyme immunoassay with human native histone H2B isolated from the eukaryotic celline HeLa for the quantitative and qualitative detection of antibodies against histone H2B in human serum.

The assay is a tool in the differential diagnosis of systemic lupus erythematosus (SLE).

## 2 Clinical Application and Principle of the Assay

Histones are small basic proteins of 11 to 21 kDa constituting a set of interacting proteins that organize and constrain the topology of DNA in most eukaryotic cells into a particle called nucleosome. In this particle the four histones H2A, H2B, H3 and H4 form an octamer (H2A, H2B, H3, H4)<sub>2</sub> around which 146 bp of dsDNA is wrapped by 1 ¾ turns. Histone H1 interacts with the nucleosome and together with linked-DNA connects neighboring nucleosomes.

Antibodies against histones occur in patients with SLE with a frequency of 50 %, during active phases of SLE with an even higher frequency (80%).

They are not specific for SLE, especially when not associated with anti-dsDNA antibodies. They are also found in other inflammatory rheumatoid diseases, like rheumatoid arthritis (low titers and mainly against H1), or systemic sclerosis and juvenile chronic arthritis. Moreover these antibodies are found in primary biliary cirrhosis (up to 70%) and autoimmune hepatitis (up to 35%), too.

Anti-histone antibodies, in the absence of antibodies to dsDNA, do play an important role in the diagnosis of drug-induced LE being found in 95% of cases. These drugs are inducing joint and skin symptoms, as well as the different anti-histone specificities, while never renal manifestations. Anti-histone antibodies in drug-induced LE are temporary and disappear within a few months after withdrawal of the inducing drug.

### Principle of the test

Serum samples diluted 1:101 are incubated in the microplates coated with the specific antigen. Patient's antibodies, if present in the specimen, bind to the antigen. The unbound fraction is washed off in the following step. Afterwards anti-human immunoglobulins conjugated to horseradish peroxidase (conjugate) are incubated and react with the antigen-antibody complex of the samples in the microplates. Unbound conjugate is washed off in the following step. Addition of TMB-substrate generates an enzymatic colorimetric (blue) reaction, which is stopped by diluted acid (color changes to yellow). The intensity of color formation from the chromogen is a function of the amount of conjugate bound to the antigen-antibody complex and this is proportional to the initial concentration of the respective antibodies in the patient sample.

### 3 Kit Contents

<b>TO BE RECONSTITUTED</b>				
Item	Quantity	Cap color	Solution color	Description / Contents
Sample Buffer (5x)	1 x 20ml	White	Yellow	5 x concentrated Tris, sodium chloride (NaCl), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Wash Buffer (50x)	1 x 20ml	White	Green	50 x concentrated Tris, NaCl, Tween 20, sodium azide < 0.1% (preservative)
<b>READY TO USE</b>				
Item	Quantity	Cap color	Solution color	Description / Contents
Negative Control	1 x 1.5ml	Green	Colorless	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Positive Control	1 x 1.5ml	Red	Yellow	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Cut-off Calibrator	1 x 1.5ml	Blue	Yellow	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Calibrators	6 x 1.5ml	White	Yellow *	Concentration of each calibrator: 0, 3, 10, 30, 100, 300 U/ml. Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Conjugate, IgG	1 x 15ml	Blue	Blue	Containing: Anti-human immunoglobulins conjugated to horseradish peroxidase, bovine serum albumin (BSA)
TMB Substrate	1 x 15ml	Black	Colorless	Stabilized tetramethylbenzidine and hydrogen peroxide (TMB/H <sub>2</sub> O <sub>2</sub> )
Stop Solution	1 x 15ml	White	Colorless	1M Hydrochloric Acid
Microtiter plate	12 x 8 well strips	N/A	N/A	With breakaway microwells. Refer to paragraph 1 for coating.
* Color increasing with concentration				
<b>MATERIALS REQUIRED, BUT NOT PROVIDED</b>				
Microtiter plate reader 450 nm reading filter and recommended 620 nm reference filter (600-690 nm). Glass ware (cylinder 100-1000ml), test tubes for dilutions. Vortex mixer, precision pipettes (10, 100, 200, 500, 1000 µl) or adjustable multipipette (100-1000µl). Microplate washing device (300 µl repeating or multichannel pipette or automated system), adsorbent paper. Our tests are designed to be used with purified water according to the definition of the United States Pharmacopeia (USP 26 - NF 21) and the European Pharmacopeia (Eur.Ph. 4th ed.).				

### 4 Storage and Shelf Life

Store all reagents and the microplate at 2-8 °C / 35.6-46.4 °F, in their original containers. Once prepared, reconstituted solutions are stable at 2-8 °C / 35.6-46.4 °F for 1 month. Reagents and the microplate shall be used within the expiry date indicated on each component, only. Avoid intense exposure of TMB solution to light. Store microplates in designated foil, including the desiccant, and seal tightly.

## 5 Precautions of Use

### 5.1 Health hazard data

**THIS PRODUCT IS FOR IN VITRO DIAGNOSTIC USE ONLY.** Thus, only staff trained and specially advised in methods of in vitro diagnostics may perform the kit. Although this product is not considered particularly toxic or dangerous in conditions of the intended use, refer to the following for maximum safety:

#### ***Recommendations and precautions***

This kit contains potentially hazardous components. Though kit reagents are not classified being irritant to eyes and skin we recommend to avoid contact with eyes and skin and wear disposable gloves.

**WARNING !** Calibrators, Controls and Buffers contain sodium azide ( $\text{NaN}_3$ ) as a preservative.  $\text{NaN}_3$  may be toxic if ingested or adsorbed by skin or eyes.  $\text{NaN}_3$  may react with lead and copper plumbing to form highly explosive metal azides. On disposal, flush with a large volume of water to prevent azide build-up. Please refer to decontamination procedures as outlined by CDC or other local/national guidelines.

**Do not smoke, eat or drink when manipulating the kit. Do not pipette by mouth.**

All human source material used for some reagents of this kit (controls, standards e.g.) has been tested by approved methods and found negative for HbsAg, Hepatitis C and HIV 1. However, no test can guarantee the absence of viral agents in such material completely. Thus handle kit controls, standards and patient samples as if capable of transmitting infectious diseases and according to national requirements.

The kit contains material of animal origin as stated in the table of contents, handle according to national requirements.

### 5.2 General directions for use

In case that the product information, including the labeling, is defective or incorrect please contact the manufacturer or the supplier of the test kit.

Do not mix or substitute Controls, Calibrators, Conjugates or microplates from different lot numbers. This may lead to variations in the results.

Allow all components to reach room temperature (20-32 °C / 68-89.6 °F) before use, mix well and follow the recommended incubation scheme for an optimum performance of the test.

**Incubation: We recommend test performance at 30 °C / 86 °F for automated systems.**

Never expose components to higher temperature than 37 °C / 98.6 °F.

Always pipette substrate solution with brand new tips only. Protect this reagent from light. Never pipette conjugate with tips used with other reagents prior.

**A definite clinical diagnosis should not be based on the results of the performed test only, but should be made by the physician after all clinical and laboratory findings have been evaluated. The diagnosis is to be verified using different diagnostic methods.**

## 6 Sample Collection, Handling and Storage

---

Use preferentially freshly collected serum samples. Blood withdrawal must follow national requirements. Do not use icteric, lipemic, hemolysed or bacterially contaminated samples. Sera with particles should be cleared by low speed centrifugation (<1000 x g). Blood samples should be collected in clean, dry and empty tubes.

After separation, the serum samples should be used during the first 8h, respectively stored tightly closed at 2-8 °C / 35.6-46.4 °F up to 48h, or frozen at - 20 °C / - 4 °F for longer periods

## 7 Assay Procedure

---

### 7.1 Preparations prior to starting

#### **Dilute concentrated reagents:**

Dilute the concentrated sample buffer 1:5 with distilled water (e.g. 20 ml plus 80 ml).

Dilute the concentrated wash buffer 1:50 with distilled water (e.g. 20 ml plus 980 ml).

To avoid mistakes we suggest to mark the cap of the different calibrators.

#### **Samples:**

Dilute serum samples 1:101 with sample buffer (1x)

e.g. 1000 µl sample buffer (1x) + 10 µl serum. Mix well !

#### **Washing:**

Prepare 20 ml of diluted wash buffer (1x) per 8 wells or 200 ml for 96 wells

e.g. 4 ml concentrate plus 196 ml distilled water.

#### **Automated washing:**

Consider excess volumes required for setting up the instrument and dead volume of robot pipette.

#### **Manual washing:**

Discard liquid from wells by inverting the plate. Knock the microwell frame with wells downside vigorously on clean adsorbent paper. Pipette 300 µl of diluted wash buffer into each well, wait for 20 seconds. Repeat the whole procedure twice again.

#### **Microplates:**

Calculate the number of wells required for the test. Remove unused wells from the frame, replace and store in the provided plastic bag, together with desiccant, seal tightly (2-8 °C / 35.6-46.4 °F)

## 7.2 Pipetting Scheme

We suggest pipetting calibrators, controls and samples as follows:

For <i>QUANTITATIVE</i> interpretation					For <i>QUALITATIVE</i> interpretation				
	1	2	3	4...		1	2	3	4...
<b>A</b>	Cal A	Cal E	P1		<b>A</b>	NC	P2		
<b>B</b>	Cal A	Cal E	P1		<b>B</b>	NC	P2		
<b>C</b>	Cal B	Cal F	P2		<b>C</b>	CC	P3		
<b>D</b>	Cal B	Cal F	P2		<b>D</b>	CC	P3		
<b>E</b>	Cal C	PC	P3		<b>E</b>	PC	...		
<b>F</b>	Cal C	PC	P3		<b>F</b>	PC	...		
<b>G</b>	Cal D	NC	...		<b>G</b>	P1	...		
<b>H</b>	Cal D	NC	...		<b>H</b>	P1	...		

CalA: calibrator A

CalB: calibrator B

CalC: calibrator C

CalD: calibrator D

CalE: calibrator E

CalF: calibrator F

PC: positive control

NC: negative control




CC: cut-off calibrator

P1: patient 1

P2: patient 2

P3: patient 3

## 7.3 Test Steps

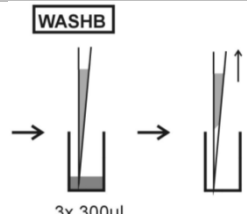Step	Description
1.	Ensure preparations from step 7.1 above have been carried out prior to pipetting.
2.	Use the following steps in accordance with quantitative/ qualitative interpretation results desired:
<b>CONTROLS &amp; SAMPLES</b>	
3.	 <p>Pipette into the designated wells as described in chapter 7.2 above, 100 µl of either:</p> <ol style="list-style-type: none"> <li>Calibrators (CAL.A to CAL.F) for <i>QUANTITATIVE</i> or</li> <li>Cut-off Calibrator (CC) for <i>QUALITATIVE</i> interp.</li> </ol> <p>and 100 µl of each of the following:</p> <ul style="list-style-type: none"> <li>Negative control (NC) and Positive control (PC), and</li> <li>Patients diluted serum (P1, P2...)</li> </ul>
4.	 <p>Incubate for 30 minutes at 20–32 °C / 68–89.6 °F.</p>
5.	 <p>Wash 3x with 300 µl washing buffer (diluted 1:50).</p>



Product Ref.	3153
Product Desc.	Histone-H2B
Manual Rev. No.	004 : 2024-12-13

### CONJUGATE


6.

CONJ



Pipette 100 µl conjugate into each well.


7.

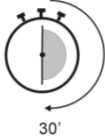

Incubate for 30 minutes at 20–32 °C / 68–89.6 °F.

8.

WASHB

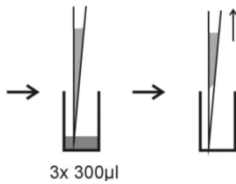

Wash 3x with 300 µl washing buffer (diluted 1:50).

### SUBSTRATE

9.

SUB



Pipette 100 µl TMB substrate into each well.

10.




Incubate for 30 minutes at 20 – 32 °C / 68 – 89.6 °F, protected from intense light.

### STOP

11.

STOP

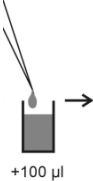

Pipette 100 µl stop solution into each well, using the same order as pipetting the substrate.

12.




Incubate 5 minutes minimum.

13.

Agitate plate carefully for 5 sec.

14.

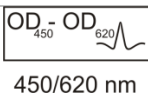

Read absorbance at 450 nm (recommended 450/620 nm) within 30 minutes.



Product Ref.	3153
Product Desc.	Histone-H2B
Manual Rev. No.	004 : 2024-12-13

## 8 Quantitative and Qualitative Interpretation

For **quantitative interpretation** establish the standard curve by plotting the optical density (OD) of each calibrator (y-axis) with respect to the corresponding concentration values in U/ml (x-axis). For best results we recommend log/lin coordinates and 4-Parameter Fit. From the OD of each sample, read the corresponding antibody concentrations expressed in U/ml.

Normal Range	Equivocal Range	Positive Results
< 12 U/ml	12 - 18 U/ml	>18 U/ml

### Example of a standard curve

**Do NOT use this example for interpreting patient's result**

Calibrators IgG	OD 450/620 nm	CV % (Variation)
0 U/ml	0.060	2.8
3 U/ml	0.137	1.0
10 U/ml	0.274	2.1
30 U/ml	0.572	2.7
100 U/ml	1.327	0.2
300 U/ml	2.289	0.2

### Example of calculation

Patient	Replicate (OD)	Mean (OD)	Result (U/ml)
P 01	0.861/0.892	0.877	53.9
P 02	0.569/0.567	0.568	28.5

Samples above the highest calibrator range should be reported as >Max. They should be diluted as appropriate and re-assayed. Samples below calibrator range should be reported as < Min.

For lot specific data, see enclosed quality control leaflet. Medical laboratories might perform an in-house quality control by using own controls and/or internal pooled sera, as foreseen by national regulations.

Each laboratory should establish its own normal range based upon its own techniques, controls, equipment and patient population according to their own established procedures.

In case that the values of the controls do not meet the criteria the test is invalid and has to be repeated.

The following technical issues should be verified: Expiration dates of (prepared) reagents, storage conditions, pipettes, devices, photometer, incubation conditions and washing methods.

If the items tested show aberrant values or any kind of deviation or that the validation criteria are not met without explicable cause please contact the manufacturer or the supplier of the test kit.

For **qualitative interpretation** read the optical density of the cut-off calibrator and the patient samples. Compare patient's OD with the OD of the cut-off calibrator. For qualitative interpretation we recommend to consider sera within a range of 20% around the cut-off value as equivocal. All samples with higher ODs are considered positive, samples with lower ODs are considered negative.

**Negative:** OD patient < 0.8 x OD cut-off  
**Equivocal:** 0.8 x OD cut-off ≤ OD patient ≤ 1.2 x OD cut-off  
**Positive:** OD patient > 1.2 x OD cut-off



Product Ref.	3153
Product Desc.	Histone-H2B
Manual Rev. No.	004 : 2024-12-13

## 9 Technical Data

Sample material:	serum
Sample volume:	10 µl of sample diluted 1:101 with 1x sample buffer
Total incubation time:	90 minutes at 20–32 °C / 68-89,6 °F
Calibration range:	0-300 U/ml
Analytical sensitivity:	1.0 U/ml
Storage:	at 2–8 °C / 35.6–46.4 °F use original vials only.
Number of determinations:	96 tests

## 10 Performance Data

### 10.1 Analytical sensitivity

Testing sample buffer 30 times on **AESKULISA® Histone-H2B** gave an analytical sensitivity of 1.0 U/ml.

### 10.2 Specificity and sensitivity

The microplate is coated with native human histone H2B. No crossreactivities to other autoantigens have been found. Antibodies against histones occur in patients with SLE with a frequency of 50 %, during active phases of SLE with an even higher frequency (80 %).

### 10.3 Linearity

Chosen sera have been tested with this kit and found to dilute linearly. However, due to the heterogeneous nature of human autoantibodies there might be samples that do not follow this rule.

Sample No.	Dilution Factor	Measured (U/ml)	Expected (U/ml)	Recovery (%)
1	1 / 100	48.8	45.0	108.4
	1 / 200	21.3	22.5	94.6
	1 / 400	10.9	11.3	96.5
	1 / 800	5.8	5.6	103.5
2	1 / 100	28.0	26.3	107.0
	1 / 200	12.4	13.2	93.9
	1 / 400	5.9	6.6	90.0
	1 / 800	3.5	3.3	106.0

## 10.4 Precision

To determine the precision of the assay, the variability (intra and inter-assay) was assessed by examining its reproducibility on three serum samples selected to represent a range over the standard curve.

Intra-assay		
Sample No.	Mean (U/ml)	CV (%)
1	54.9	2.2
2	31.2	2.1
3	10.5	0.7

Inter-assay		
Sample No.	Mean (U/ml)	CV (%)
1	51.8	6.6
2	26.5	4.0
3	8.0	0.9

## 10.5 Calibration

Due to the lack of international reference calibration this assay is calibrated in arbitrary units (U/ml).

## 11 Literature

**Gioud M, Ait Kaci A, Monier JC (1982).** Histone antibodies in systemic lupus erythematosus. A possible diagnostic tool. *Arthritis Rheum* 25: 407-413.





















**Rubin RL, Joslin FG, Tan EM (1982).** Specificity of anti-histone antibodies in systemic lupus erythematosus. *Arthritis Rheum* 25: 779-782.

**Shoenfeld Y and Segol O (1989).** Anti-histone antibodies in SLE and other autoimmune diseases. *Clin and Exp Rheum* 7: 265-271.

**Cohen MG and Welb J (1989).** Anti-histone antibodies in rheumatoid arthritis and Felty's syndrome. *Arthr. Rheum* 32: 10

**Burlingame RW (1997).** The clinical utility of anti-histone antibodies. *Clinics in Laboratory Medicine* 17: 367-378.

## 12 Regulatory Symbols

	- Diagnosi in vitro - Pour diagnostic in vitro - In Vitro Diagnostikum - Para uso Diagnóstico in vitro	- For in vitro diagnostic use - Para uso diagnóstico in vitro - In Vitro Διαγνωστικό μέσο
	- Numero d'ordine - Référence Catalogue - Bestellnummer - Número de catálogo	- Catalogue number - Numéro de catálogo - Αριθμός παραγγελίας
	- Descrizione lotto - Lot - Chargen Bezeichnung - Lote	- Lot - Lote - Χαρακτηρισμός παρτίδας
	- Identificatore univoco del dispositivo - Identifiant unique de l'appareil - eindeutige Produktidentifizierung - Identificador único do dispositivo	- Unique Device Identifier - Identificador único del dispositivo - Μοναδικό αναγνωριστικό συσκευής
	- Conformità europea - Déclaration CE de Conformité - Europäische Konformität - Declaração CE de Conformidade	- EC Declaration of Conformity - Declaración CE de Conformidad - Ευρωπαϊκή συμφωνία
	- 96 determinazioni - 96 tests - 96 Bestimmungen - 96 Testes	- 96 tests - 96 pruebas - 96 προσδιορισμοί
	- Respectare le istruzioni elettroniche per l'uso - Voir les instructions d'utilisation électronique - Elektronische Gebrauchsanweisung beachten - Seguir as instruções electrónicas de utilização	- See electrical instructions for use - Siga las instrucciones electrónicas de uso - Ακολουθήστε τις ηλεκτρονικές οδηγίες χρήσης
	- Da utilizzarsi entro - Utilise avant le - Verwendbar bis - Utilizar antes de	- Use by - Utilizar antes de - Χρήση μέχρι
	- Conservare a 2-8°C (35.6-46.4°F) - Conserver à 2-8°C (35.6-46.4°F) - Lagerung bei 2-8°C (35.6-46.4°F) - Conservar entre 2-8°C (35.6-46.4°F)	- Store at 2-8°C (35.6-46.4°F) - Conservar a 2-8°C (35.6-46.4°F) - Φυλάσσεται στους 2-8°C (35.6-46.4°F)
	- Prodotto da - Fabriqué par - Hergestellt von - Fabricado por	- Manufactured by - Fabricado por - Κατασκευάζεται από
	- Calibratore cut-off - Etalon Seuil - Grenzwert Kalibrator - Calibrador de cut-off	- Cut off Calibrator - Calibrador de cut-off - Οριακός ορός Αντιδραστήριο βαθμονόμησης
	- Controllo positivo - Contrôle Positif - Positiv Kontrolle - Controllo positivo	- Positive Control - Control Positivo - Θετικός ορός ελέγχου
	- Controllo negativo - Contrôle Négatif - Negativ Kontrolle - Controllo negativo	- Negative Control - Control Negativo - Αρνητικός ορός ελέγχου
	- Calibratore - Etalon - Kalibrator - Calibrador	- Calibrator - Calibrador - Αντιδραστήριο βαθμονόμησης
	- Coniugato - Conjugé - Konjugat - Conjugado	- Conjugate - Conjugado - Σύζευγμα
	- Micropiastra rivestita - Microplaque sensibilisée - Beschichtete Mikrotiterplatte - Microplaca revestida	- Coated microtiter plate - Microplaca sensibilizada - Επικαλυμμένη μικροπλάκα
	- Tampone di lavaggio - Tampon de Lavage - Waschpuffer - Solução de lavagem	- Wash buffer - Solución de lavado - Ρυθμιστικό διάλυμα πλύσης
	- Tampone substrato - Substrat - Substratpuffer - Substrato	- Substrate buffer - Tampón sustrato - Ρυθμιστικό διάλυμα υποστρώματος
	- Reagente bloccante - Solution d'Arrêt - Stopreagenz - Solução de paragem	- Stop solution - Solución de parada - Αντιδραστήριο διακοπής αντίδρασης
	- Tampone campione - Tampon Echantillons - Probenpuffer - Diluente de amostra	- Sample buffer - Tampón Muestras - Ρυθμιστικό διάλυμα δειγμάτων